

Controlling nature just got easier



LBBTM **Lombardi**TM
Bionator Bite Appliance

MASEL[®]

LBB™ Lombardi™

Bionator Bite Appliance

The Lombardi Bionator Bite Appliance Delivers

The *LBB*™ delivers what many retainers and custom appliances deliver, without the impression, wait, or cost. It will open the bite in deep bite cases, aid in the treatment of TMJ disorders and bruxism, and facilitate tooth movement during orthodontic treatment.

Include One in Your Deep Bite Patient Treatment Kit

The *LBB* can assist your treatment method, for nearly all patients, with varying degrees of malocclusion. Plus the appliance comes in one size, fitting patients from 8 to 80 years of age. Each *Lombardi*™ Bionator Bite Appliance comes with a storage case to keep the device flat when not being worn. This helps to “recharge” the spring-action hinges. Include the *LBB* as part of your deep bite treatment kit.

Patient Comfort Equals Patient Compliance

The most dramatic bite improvement occurs when a patient wears the appliance all night while asleep and for 3 to 4 hours during quiet daytime hours. Suggestions for wearing while awake include during workouts, while reading, listening to music or watching television, while online, or doing homework. It is simple to clean with soap and water, or with a toothbrush and toothpaste.

“It’s saving time and reducing costs; something patients really appreciate.”

Dr. Steven Nederlow, Board Certified Orthodontist in private practice in Minnesota

Features:

- Ready to go, no impressions needed
- No risk of aspiration
- Compatible with all systems, brackets, and techniques
- Assists the function of elastics, cross elastics, and headgear
- Reduces bond failures
- Opens the bite in deep-bite cases and protects occlusion
- Aids in the treatment of TMJ disorders
- Cost effective and highly profitable
- Latex-free

Orthodontic treatment can progress more effectively when you include the LBB Appliance.

Drew R. Lombardi, D.M.D. invented the *LBB* Appliance to minimize the negative effects of bruxing, clenching, and grinding that interfere with orthodontic treatment, so tooth movement proceeds efficiently.



Drew R. Lombardi, D.M.D.

Dr. Drew Lombardi is the inventor of the *Lombardi* Bionator Bite Appliance.


Dr. Lombardi has two private orthodontic practices in Fort Lee and Jersey City, NJ. A Diplomate of the American Board of Orthodontics, he has taught and lectured at the Jersey City Medical Center, and Fairleigh Dickinson University, School of Dentistry, Department of Orthodontics, where he was a Clinical Assistant Professor in Orthodontics for years until the school's closing.

He was selected by a poll of over 4,500 New Jersey dentists as the best orthodontist in Bergen and Hudson Counties, and among the top 10 in the entire state of New Jersey, as featured in New Jersey Monthly Magazine, July 2002. He was again selected to this list in 2005, 2009, and 2010. (Each year that it's been done!)

www.LBB1.com

Bite pressure is transferred from posterior to anterior region

- More sensitive anteriors signal brain to decrease clenching
- Helps open bite

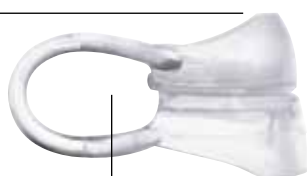

Molar region is left free of occlusal pressure

- Allows posterior teeth to extrude naturally when worn regularly

Flexible channel design fits loosely over anterior sextants

- Device can be worn comfortably by most patients
- Allows for malocclusion and overjet
- Accommodates orthodontic brackets, wires, elastics, tooth movement, and case progress without adjustment
- Cut out for maxillary frenum to ensure comfort




Flexible hinge applies gentle pressure in opposite, bite-opening directions

- Holds the appliance in place so mandible rests comfortably, not where the malocclusion normally forces it to go

Anterior occlusal plane provides protection for upper and lower anterior teeth

- Keeps lower anteriors from hitting lingual surface of uppers
- Protects upper incisal edges from hitting lower brackets

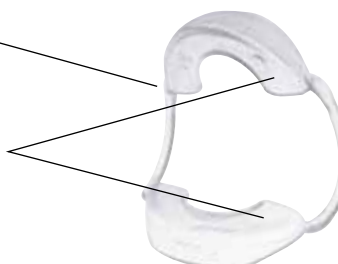

"It's a fast, easy, inexpensive way to help patients in a number of situations I encounter every day."

Dr. Howard S. Glazer—General Practitioner, practicing in New Jersey

"Even non-compliant patients will wear the LBB. We're seeing much faster orthodontic treatment ... and happier patients."

Dr. Richard Isaacson—Orthodontist, in family practice in New Jersey



Kit Includes:

- 1 LBB
- 1 Protective Case
- Patient Care Instructions

DESCRIPTION	ITEM NO.
LBB Appliance Kit	5200-100



LBB™ Lombardi™ *Bionator Bite Appliance*

MASEL®

Ortho Organizers, Inc.
1822 Aston Avenue
Carlsbad, CA 92008
U.S.A.
MaselOrtho.com

To Order by Phone

Domestic

Direct:

Ph 760 448 8600
Fax 760 448 8602

Toll-free:

Ph 888 851 0533
Fax 800 888 7244

Email Sales@Masel.net

International

Direct:

Ph +(1) 215 458 9140 +(1) 760 448 8600
Fax +(1) 215 458 9155 +(1) 760 448 8607

Toll-free:

Australia

Ph 0011 800 0986 2735
Fax 0011 800 0996 2735

Austria, Belgium, Denmark, Germany, Italy, The Netherlands, Norway, Switzerland, UK

Ph 00 800 0986 2735
Fax 00 800 0996 2735

Canada

Ph +(1) 866 510 4252
Fax +(1) 866 510 1822

Mexico

Ph 001 866 510 4252
Fax 001 866 510 1822

South Africa

Ph 09 800 0986 2735
Fax 09 800 0996 2735

To Order Online

Order.MaselOrtho.com