

CV SuperForce™ Archwires

Among the strongest tooth-moving forces for thermal nickel titanium archwires

At room temperature, CV SuperForce thermal nickel titanium wire is consistently soft, and makes ligation easier for you and more comfortable for the patient. Highly responsive to chilling, the wire can be softened even further so the most difficult cases can be ligated quickly and easily. The force of the wire is directly between that of Elastinol™ nickel titanium SE and Therm-A-Form ULTRA™ thermal nickel titanium archwires.


Characteristics of CV SuperForce include:

- Ideal for adult patients with significant crowding, and for early- to mid-stage treatment with moderate to severe crowding
- Easy to engage because of the thermal properties
- Light and continuous forces, utilizing shape memory and super elastic property to move teeth
- Etched center marks, three marks on uppers and one mark on lowers
- Heat activation transforms the wire to a super elastic state at 65-75° F (18-24° C)
- Available in the Natural archform
- 10 per pack




Sizes		Natural	
INCHES	MM	UPPER	LOWER
.014	.36	4805-114	4805-014
.016	.40	4805-116	4805-016
.018	.46	4805-118	4805-018
.016 x .016	.40 x .40	4805-161	4805-061
.016 x .022	.40 x .56	4805-162	4805-062
.017 x .025	.43 x .64	4805-172	4805-072
.018 x .025	.46 x .64	4805-182	4805-082
.019 x .025	.48 x .64	4805-190	4805-090
.021 x .025	.53 x .64	4805-121	4805-021

Ultimate Wireforms, Inc.
Nickel Titanium Wire Force Chart
UltraTherm Plus
Round Wires

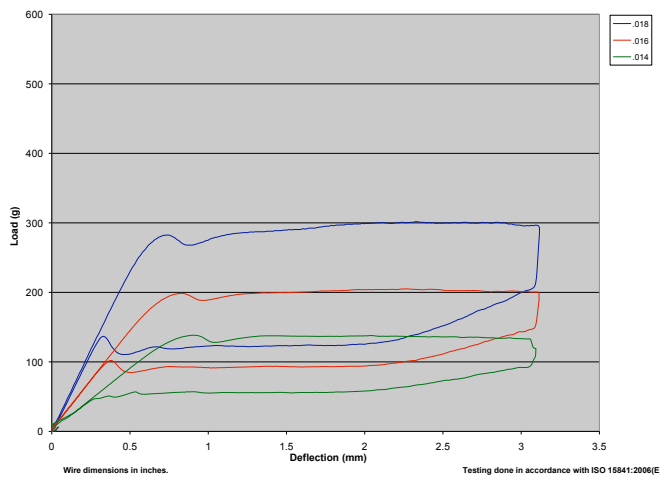

Ultimate Wireforms, Inc.
Nickel Titanium Wire Force Chart
UltraTherm Plus
Rectangular & Square Wires

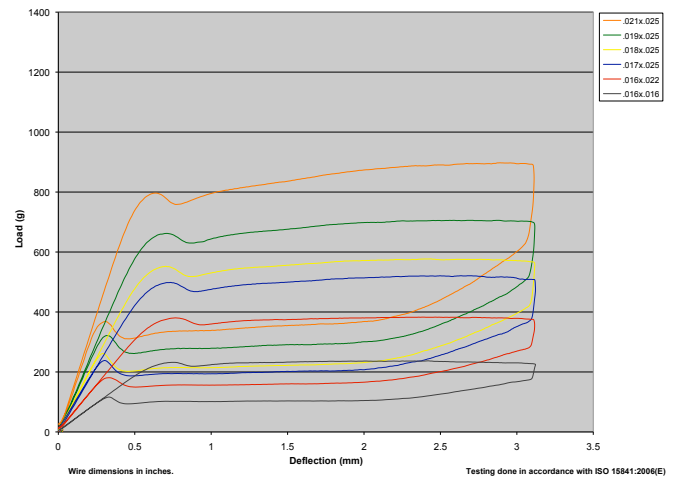

To order by phone

In the U.S.

Toll-free: 888 851 0533
Fax: 800 888 7244

Outside the U.S.

Direct: +(1) 760 448 8600
Fax: +(1) 760 448 8607

Visit our Website at MaselOrtho.com