


**ADHESIVES
AND CEMENTS**



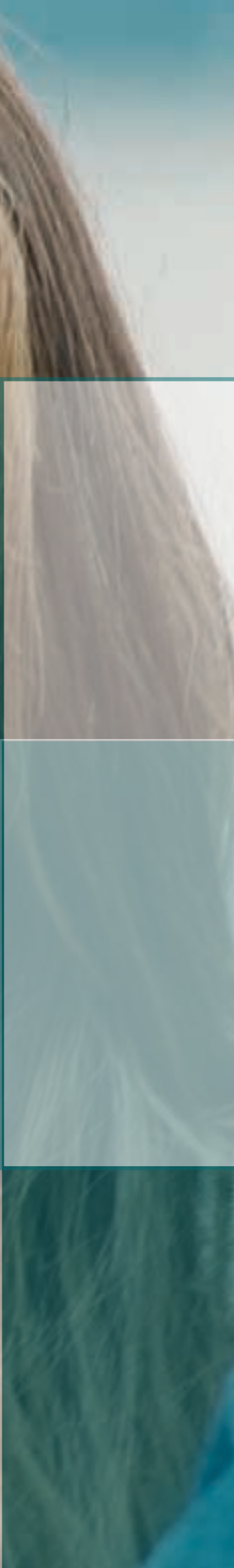

Everest™ Adhesives **D-3**
Everest Cements and Bracket Conditioner **D-6**

Everest Light-Cure Orthodontic Adhesive

The smooth, tacky composite resin delivers ideal handling characteristics which eliminates bracket float and allows for easy flash clean up. Adhesive provides the optimum properties for bonding metal, ceramic, or plastic brackets using light-cure technology. Light-activation provides accurate bracket placement ensured by allowing for sufficient working time. Produces a strong light-cured chemical bond which also debonds predictably. Archwires can be placed immediately after fully curing the brackets. Contains fluoride to help protect against potential decalcification.

Everest Light-Cure Kit in Syringes*

Kit includes:

- 4 Push Syringes (3.5 grams ea.)
- 1 Gel Etch Syringe (4 grams)
- 1 Sealant Resin (7 ml)
- 5 Syringe Tips (22-Gauge)
- 2 Mixing Wells (10-Cavity)
- 1 Mixing Pad
- 50 Brushes
- 1 Instruction For Use Sheet

ITEM NO.
490-300



Everest Light-Cure Kit in Capsules*

Kit includes:

- 40 Capsules (.35 grams each)
- 1 Sealant Resin (7 ml)
- 1 Green Liquid Etchant (9 grams)
- 1 Mixing Pad
- 50 Brushes
- 1 Instruction For Use Sheet



ITEM NO.
490-320

Everest Dispensing Gun

1 Capsule Gun



ITEM NO.
490-220

Everest Light-Cure Kit*

Kit includes:

- 1 Push Syringe (3.5 grams)
- 1 Sealant Resin (3 ml)
- 1 Green Gel Etchant Syringe (4 grams)
- 5 Syringe Tips (22-Gauge)
- 5 Bonding Brushes
- 1 Instruction For Use Sheet

ITEM NO.
490-310



Everest Light-Cure Reorder Components

DESCRIPTION	ITEM NO.
4 Light-Cure Adhesive Syringes Pack (3.5 grams ea.)	490-111
1 Light-Cure Adhesive Single Syringe (3.5 grams ea.)	490-112
40 Light-Cure Single Use Capsules (.35 grams ea.)	490-221
1 Sealant Resin (7 ml)	490-200
1 Gel Etchant Large Bottle (330 grams)*	490-209
1 Gel Etchant Small Bottle (57 grams)*	490-210
5 Mixing Pads	490-505
1 Brush Holder	490-504
100 Disposable Brushes	490-510
50 Disposable Long-Handled Brushes	490-520



* Ground shipment recommended – This product is classified as hazardous material as per the U.S. Department of Transportation Pipeline and Hazardous Materials Safety Administration. Extra shipping charges may apply.

EVEREST
Orthodontic Adhesive

Everest Low Viscosity, Light-Cure Flowable Adhesive*

Everest Low Viscosity Light-Cure Flowable Adhesive features a precision-controlled syringe which prevents overflow and run-on. It ensures accuracy as it flows easily into difficult areas, but is reliable not to slump or run until light curing occurs. The disposable tip provides precise application, accessibility, and control, resulting in excellent aesthetics. It is dependable, strong, and more resistant to shock and stress. This adhesive is ideal for bracket bases with greater than an 80-gauge mesh and can be used in many ways, such as bonding lingual retainers, Invisalign® attachments, occlusal build-ups for opening a posterior bite, indirect bonding, and much more! Contains fluoride to help protect against potential decalcification.

Everest Low Viscosity Adhesive Kit*

Kit Includes:

- 3 Flowable Adhesive Syringes (2 grams ea.)
- 1 Green Etchant Syringe (1.6 grams)
- 1 Light-Cure Enamel Bonding Sealant Resin (3 ml)
- 25 Flowthru 22 Gauge Syringe Tips
- 10 Bonding Brushes
- 1 Instruction For Use Sheet

ITEM NO.
490-330

Everest Low Viscosity Reorder Components

DESCRIPTION	ITEM NO.
1 Low Viscosity Single Syringe (2 grams) with 10 Flowthru 22-Gauge Syringe Tips	490-223
20 Low Viscosity Single Use Capsules (.35 grams ea.)	490-222
1 Dispensing Gun Capsule Gun	490-220
25 Syringe Tips (Gel Etch/Low Viscosity Adhesive)	490-235
1 Sealant Resin (7 ml)	490-200
1 Gel Etchant Large Bottle (330 grams)*	490-209
1 Gel Etchant Small Bottle (57 grams)*	490-210
1 Brush Holder	490-504
100 Disposable Brushes	490-510
5 Disposable Long-Handled Brushes	490-520



* Ground shipment recommended – This product is classified as hazardous material as per the U.S. Department of Transportation Pipeline and Hazardous Materials Safety Administration. Extra shipping charges may apply.

*Call for availability.

EVEREST
Orthodontic Adhesive

Everest No-Mix Orthodontic Adhesive (Primer Activated) Bonding System

Provides the optimum properties for bonding ceramic, metal, or plastic brackets using a primer activated no-mix bonding system. This is a chemically-cured system which does not require pre-mixing. The thick, tacky adhesive resin delivers ideal handling characteristics which eliminate bracket flotation and allow for easy flash clean up. Archwires can be placed 5 minutes after bonding the last bracket. Contains fluoride to help protect against potential decalcification.

Everest No-Mix Orthodontic Adhesive Kit in Syringes (Primer-Activated)*

Kit Includes:

- 4 Push Syringes (3.5 grams ea.)
- 1 Enamel Catalyst Primer (9 ml)
- 1 Gel Etchant Syringe (4 grams)
- 5 Syringe Tips (22-Gauge)
- 2 Mixing Wells (10-Cavity)
- 1 Mixing Pad
- 50 Brushes
- 1 Instruction For Use Sheet



ITEM NO.
490-340

Everest No-Mix Orthodontic Adhesive Mini Kit in Syringe (Primer-Activated)*

Kit Includes:

- 1 Push Syringe (3.5 grams)
- 1 Enamel Catalyst Primer (3 ml)
- 1 Gel Etchant Syringe (4 grams)
- 5 Syringe Tips (22-Gauge)
- 5 Bonding Brushes
- 1 Instruction For Use Sheet



ITEM NO.
490-350

Everest Etchant Bonding Supplies

Gel stays where placed without overflowing and provides highly visible contrast against tooth surfaces. Non-drying water soluble consistency with ideal medium viscosity. Traditional 37% Phosphoric acid.

Everest Gel Etch Kit, Economy Size*

Kit Includes:

- 4 Syringes (4 grams ea.)
- 20 Plastic Flowthru Tips (22-Gauge)



ITEM NO.
490-360

Everest Gel Etch Mini Kit*

Kit Includes:

- 1 Syringe (4 grams)
- 5 Plastic Flowthru Tips (22-Gauge)



ITEM NO.
490-370

Everest Etchant Reorder Components

DESCRIPTION	ITEM NO.
1 Gel Etchant Large Bottle (330 grams)*	490-209
1 Gel Etchant Small Bottle (57 grams)*	490-210
1 Liquid Gel Etchant Large Bottle (29.6 ml)*	490-211
1 Liquid Gel Etchant Small Bottle (7 ml)*	490-212
25 Syringe Tips (Gel Etch/Low Viscosity Adhesive)	490-235



* Ground shipment recommended – This product is classified as hazardous material as per the U.S. Department of Transportation Pipeline and Hazardous Materials Safety Administration. Extra shipping charges may apply.

EVEREST
Orthodontic Adhesive

Everest Light-Cure Band Cement

Provides a longer working time and a rapid setting time with its light-activated band cement in a single paste system. Can be placed directly into bands. Excess flash is easy to clean up prior to the light curing and the cement is easy to see due to the blue color. Contains fluoride to help protect against potential decalcification.

Everest Light-Cure Orthodontic Band Cement Kit

Kit Includes:

- 6 Push Syringes (5 grams ea.)

ITEM NO.
490-390



Everest Powder/Liquid Band Cement

Everest chemically bonds to enamel without the use of etchants. Using a two-part glass ionomer band cement system of powder and liquid, it has a low film thickness and will not wash out under bands. The final setting time is 6 to 7 minutes. Contains fluoride to help protect against potential decalcification.

Everest Glass Ionomer Cement Economy Kit

Kit Includes:

- 1 Powder (100 grams)
- 2 Liquid (50 grams ea.)
- 1 Mixing Pad
- 1 Measuring Scoop
- 1 Instruction For Use Sheet



ITEM NO.
490-380

Everest Glass Ionomer Cement Starter Kit

Kit Includes:

- 1 Powder (50 grams)
- 1 Liquid (50 grams)
- 1 Mixing Pad
- 1 Measuring Scoop
- 1 Instruction For Use Sheet



ITEM NO.
490-385

Plastic Bracket Conditioner

Prepares plastic brackets for direct bonding. Brush on base prior to applying adhesive. Provides a strong bond between the bracket base and adhesive. Must be refrigerated when not in use. (10 cc)

ITEM NO.
470-170



Everest Glass Ionomer Cement Reorder Components

DESCRIPTION	ITEM NO.
1 Glass Ionomer Powder (50 grams)	490-270
1 Glass Ionomer Liquid (50 ml)	490-271

Everest Miscellaneous Reorder Components

DESCRIPTION	ITEM NO.
50 Disposable Spatulas	490-500
5 Brush Holders	490-504
5 Mixing Pads	490-505
5 Measuring Scoops	490-506
100 Disposable Brushes	490-510
50 Disposable Long-Handled Brushes	490-520
10 Mixing Wells (Dual Cavity)	490-535
1 Dispensing Gun	490-220
25 Syringe Tips (Gel Etch/Low Viscosity Adhesive)	490-235
Masel Bond Prep Primer and Conditioner	4000-803

